

CORDLESS ENDO MOTOR

CanalPro X-Move

Full Freedom of Movement





CanalPro X-Move

The wireless CanalPro X-Move endomotor provides two movement modes, a full rotating and a reciprocation motion for all common file systems. The motor stands for simple handling and great flexibility.

Features & Benefits


- Continuous rotation and a reciprocation motion for all common file systems
- Wireless endomotor for easy and flexible operation
- Integrated apex locator for safe and efficient treatments
- Insulated contra angle with super mini-head and slim neck provides convenient usage and free field of view
- Intuitive control display for easy handling
- Pre-programmed settings for COLTENE file systems saves preparation time



Normal contra angle head & neck




CanalPro X-Move contra angle super mini-head & slim neck

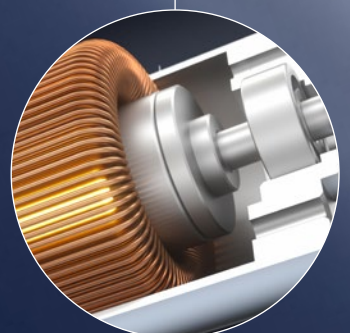

Technical Specifications

- Contra angle: 360 degrees rotation
- Torque range: 0.4 – 5.0 Ncm
- Speed range: 100 – 2500 rpm
- Adjustable reciprocating angle: 10° interval
- Adjustment range 20° – 400°
- Mini insulated contra angle with diameter of 8 mm and a height of 9.7 mm

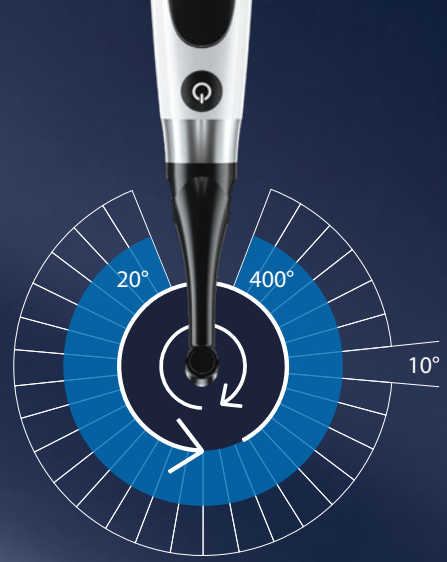

User Menu

- Pre-programmed settings for HyFlex EDM, HyFlex EDM OGSF Sequence, HyFlex CM, MicroMega One Flare, MicroMega One G, MicroMega One Curve mini, MicroMega One RECI, MicroMega 2Shape mini, Remover for HyFlex and MicroMega systems
- Included option to manually set parameters (RPM, torque, angles) for other files systems



The Motor – stronger power, lower noise

- Adjustable reciprocation angles for higher preparation efficiency
- Lower noise levels due to smooth rotation provides higher treatment comfort
- High performance brushless motor
- Control display for easy handling

Technical and settings
data of all COLTENE files:



Order Information

CanalPro X-Move Box

65002739	1 × Motor Handpiece
	1 × Charging Base
	1 × Contra Angle
	1 × Spray Nozzle
	1 × Measuring Wire
	4 × File Holder
	2 × Lip Clip
	2 × Touch Probe
	2 × "O" ring
	1 × Power Adapter
	1 × Battery

Accessories

65002742	1 × Motor Handpiece
65002743	1 × Charging Base
65002744	1 × Contra Angle
65002745	1 × Spray Nozzle
65002746	1 × Power Adapter
65002747	1 × Measuring Wire
65002748	4 × File Holder
65002749	2 × Lip Clip
65002750	2 × Touch Probe
65002751	2 × "O" ring
65002752	1 × Battery



Your Endo Guide

Practical. Innovative. Leading the way.

COLTENE provides simple, safe, and efficient solutions for every step of the endo process, carefully guiding the dentist through every endodontic situation with dedicated trainings and experienced endo consultants.

Our portfolio, which has been refined over decades, includes innovative consumables and sophisticated equipment that reliably cover all five treatment steps.



Infection
Control



Cleaning &
Shaping



Irrigation &
Drying



Obturation



Restoration

© COLTENE – www.coltene.com

Imprint / Publisher

Coltene/Whaledent AG
Feldwiesenstrasse 20
CH-9450 Altstätten
+41 71 757 5300
info.ch@coltene.com

Sales Contact

Coltene/Whaledent Ltd.
The President Suite
Kendal House, Victoria Way
Burgess Hill, West Sussex
RH15 9NF/UK
+44 1444 235486
info.uk@coltene.com

Service Center

service@coltene.com
00 800 11 0000 12 (U/FN)
Mon.-Fri. 8:00 to 17:00 MESZ

 **COLTENE**